

Multi-Cassette Autoclave



Easy:

Choose cycle, and start it, so easy!

Reliable:

Low failure rate based over 15 years' experience and customer feedback.

Gentle:

Steam is injected by positive pressured pulses, safe and gentle on expensive instruments.

Safe:

Safety valve guarantees operation.

Small:

Compact in structure, can be placed up and down the table.

Smart:

Data can be saved and traced back!

Efficient :


shortest unwrapped cycle: 6~8 minutes,
wrapped cycle: 10~15 minutes.
All Cassettes sterilized together or independently.



Save Your Time & Space!

Working Principle:

Positive pressured pulse to remove cold air inside cassette.

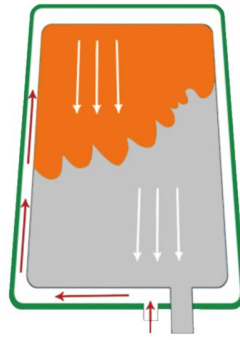

Cycle parameters

Symbol	Cycle	Sterilize Temp.	Sterilize Time	Dry Time	Positive Pulse
	Rapid (Unwrapped N)	134°C	240s	900s	1
	General (Wrapped S)	134°C	240s	3600s	3
	Porous (Unwrapped S)	134°C	600s	900s	3
	Low temp(Wrapped S)	121°C	1800s	900s	3
	Self-define (Mixed)	134°C	900s	900s	6
	Air Dry (Auxiliary)	-----	-----	900s	-----

Model Specification

Model Specifications	Celer™ 2000			Celer™ 5000			Celer™ 6000		
	S1	S2	S3	S1	S2	S3	S1	S2	S3
Length(mm)	615								
Width(mm)	320								
Height(mm)	320	470	620	320	470	620	320	470	620
Net weight(kg)	27	41	55	33	54	75	33	54	75
Power(kw)	1.3	2.6	3.9	1.3	2.6	3.9	1.3	2.6	3.9
Cassette Overall Size(L×W×Hmm)	390×196×40mm			480×196×80mm			570×196×80mm		
Cassette volume	1.8L			5.2L			6.0L		
Working Temperature	115°C ~ 134°C								
Power	AC 220V±22V; 50Hz±1Hz;								

AIKSMED

JIANGSU AIKSMED EQUIPMENT CO., LTD

ADD: Xinlian Advanced Manufacturing Industrial Park,
Daxin Town, Zhangjiagang City, Jiangsu Province, China

Email: sales@aiksmmed.com

Tel/Fax: 0086-512-5876-3098

www.steam-autoclave.com

www.aiksmmed.com



Authorized Distributor

